

# Sonosite PX

Advanced image clarity and confidence



Any patient. Anywhere. Anytime.

# See more. Do more.

## Sonosite PX

### Advanced image clarity **and** confidence

Inspired, informed, and improved by thousands of hours working side-by-side with clinicians on the front lines of medicine, Sonosite PX point-of-care ultrasound system incorporates our most advanced image clarity to give clinicians an advanced level of confidence for diagnostic and procedural applications.

**This isn't just our most advanced point-of-care ultrasound solution ever, it's YOURS.**



## Designed for the human touch

The human touch is crucial when it comes to point-of-care ultrasound. That's why Sonosite PX has been designed around every detail for the way clinicians work and deliver care.

### Touch and physical controls

Intuitive touch panel and physical buttons for quick access to the most-used controls. Unique raised buttons allow clinicians to make adjustments via feel. Sonosite PX can be operated with gloves and a sterile drape.

### Workflow efficiency features

The streamlined user interface improves workflow, allowing clinicians to find controls quickly with fewer touches. Exam types are optimized with most-needed annotations and calculations.



## Get closer at the point of care

We listened and listened some more, working with clinicians across the globe. Sonosite PX was created and inspired by the way they need to deliver care, wherever they deliver it.

### Adaptable positioning

Adaptable horizontal to vertical work surface allows clinicians to adjust Sonosite PX for optimal bedside ergonomics and bring it closer to the bedside.

### Closer access

The low-profile wheelbase fits under most beds for closer access, especially important while operating the system across the bed. Front and side levers simplify height adjustments.

#### Adjustable height

System height can be raised or lowered for a wide range of users.

#### Adaptable positioning

Work surface adapts from horizontal to vertical with a push or pull of the handle.



# Environment inspired.

## Optimized for the way you move. The way you think.

The agile, small footprint Sonosite PX stand offers maneuverability and dedicated features to keep your system organized and stocked with the essentials.

### Data management and connectivity

Wifi, Ethernet port, 1080p digital video output compatible with 1080p displays, (4) USB 2.0 ports, (2) USB 3.0 ports, (1) dedicated USB port for optional printer.<sup>†</sup>

### Durable transducer holders

Integrated holders (no more cups) for simplified clean up and ease of use.

### Triple transducer connect

Three transducer ports so your most-used transducers are at the ready with the touch of a button.

### Versatile storage solutions

The flexible storage system adapts to left or right side needs—gel holder, wipes holder, storage bin, and lockable drawer to secure supplies.

### Cord management

Power cords are confined within the column to help prevent tangles.

### Cable management system

Cable catcher and hooks help prevent tangled transducer cables. Even at the lowest stand height, cables are less likely to drag on the floor.

### Low-profile wheeled base

The low-profile wheel base helps position the system close to the bedside.



### Peripheral accessories

#### Barcode scanner

Save time by automatically linking patient data to the ultrasound scan. English only, where available.

#### Barcode expressions

Auto-populate physician ID and patient information with customized software for your Sonosite Barcode Scanner. English only, where available.

#### Sonosite ECG module

The 3-lead Sonosite ECG Module connects directly to Sonosite PX ultrasound system through a dedicated ECG port. It allows the system to display an ECG trace on the clinical monitor. Includes ECG lead wire set.

#### Medical grade printer

Add the optional Black and White Printer to print ultrasound studies.

#### Needle guide starter kits

Improve positioning to help guide needles to an exact target on the first attempt.

#### PowerPark cordless charging station


Dock the system in a central location while the battery charges for cord-free use. Cord-free use eliminates power cords dragging on the floor, risking contamination. No power cords to clean and disinfect.

<sup>†</sup> Off the stand, Sonosite PX has (2) USB 3.0 ports.




### Integrate medical informatics

Synapse Synchronicity workflow manager supports seamless credentialing, quality assurance, and billing. Compile data from all connected point-of-care ultrasound systems into a central location.

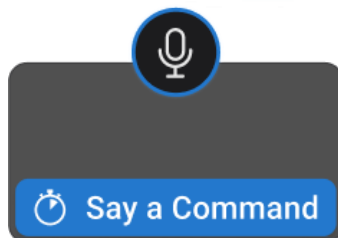


### Voice assist

Hands free system controls for common modes and features providing you the ability to complete procedures without assistance or in a sterile environment.

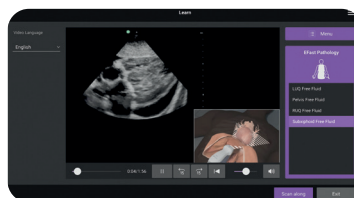

### Education at your fingertips

Sonosite is dedicated to clinician education. More than 150 scan-along education tutorials are onboard for COVID-19, diagnosis, procedures, pathology, and system and transducer use.

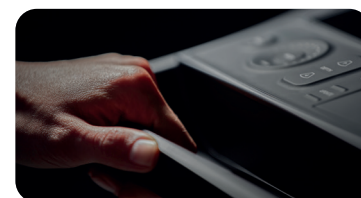

### Rest assured with data security

Sonosite PX includes robust security features to keep patient information safe and networks secure. The system meets FIPS199, FIPS200, and NIST 800-53.



# Scan with confidence

Transducers designed and manufactured by Sonosite undergo rigorous testing in-house. Our transducers aren't only durable, they're also made for superior imaging in the most demanding environments.

- IPX7 rated for immersion\*
- Drop-tested to 1 meter
- Sonosite-manufactured transducers are covered by a Sonosite 5-year standard warranty

## Optional needle guide starter kit

- CIVCO Infiniti Plus™ In-Plane Ultrasound Needle Guide
- CIVCO AccuSITE™ Out-of-Plane Ultrasound Needle Guide
- CIVCO Disposable Endocavity Needle Guide

## Available Transducers



### C5-1 ●

5-1 MHz Curved

Scan depth: 4.7–30 cm

#### Exam types:

Abdomen, Procedural Abdomen, Gynecology, Lung, Musculoskeletal, Nerve, Early Obstetrics, Obstetrics, Spine, VExUS



### C10-3 ●

10-3 MHz Curved

Scan depth: 3–18 cm

#### Exam types:

Abdomen, Gynecology, Lung, Musculoskeletal, Nerve, Spine, VExUS



### P5-1

5-1 MHz Phased

Scan depth: 5–35 cm

#### Exam types:

Abdomen, Pediatric Abdomen, Cardiac, Pediatric Cardiac, Focused Cardiac, Lung, Obstetrics, Orbital, Transcranial Doppler (TCD), VExUS



### P11-3

11-3 MHz Phased

Scan depth: 2–15 cm

#### Exam types:

Neonatal and Pediatric Cardiac, Neonatal and Pediatric Lung, Neonatal and Pediatric Abdominal, Neonatal Cephalic



### L12-3 ●

12-3 MHz Linear

Scan depth: 1.5–9 cm

#### Exam types:

Arterial, Breast, Carotid, Lung, Musculoskeletal, Nerve, Ophthalmic, PIV, Superficial, Venous



### L13-6

13-6 MHz Multi

Scan depth: 1–6 cm

#### Exam types:

Musculoskeletal, PIV, Superficial, Venous



### L15-4 ●

15-4 MHz Linear

Scan depth: 1.5–6 cm

#### Exam types:

Arterial, Breast, Carotid, Musculoskeletal, Nerve, Superficial, Venous



### L19-5 ●●

19-5 MHz Linear

Scan depth: 1–6 cm

#### Exam types:

Arterial, Lung, Musculoskeletal, Nerve, Ophthalmic, PIV, Superficial, Vascular Access, Venous



### IC10-3 ●

10-3 MHz Curved

Scan depth: 3–15 cm

#### Exam types:

Gynecology, Early Obstetrics, Obstetrics, Prostate



### T8-3\*\*

8-3 MHz Multi

Scan depth: 4–18 cm

#### Exam types:

Cardiac, Cardiac Resuscitation, Lung

\* Protected against water immersion for 30 minutes at a depth of up to 1 meter.

\*\* 1-year warranty

## Anesthesia

### Procedural Support


For clinicians working at the leading edge of ultrasound guided anesthesia and diagnostics, Sonosite PX incorporates our most advanced image clarity to help clinicians diagnose and more accurately visualize sonographic landmarks for anatomical certainty during procedures.



## Critical Care

### Caring for the Critically Ill

Allows clinicians to effectively assess cardiac function and fluid volume status with a simplified VExUS scoring exam. Sonosite PX provides cardiac functionality that supports all levels of users, from novice to advanced. The streamlined user interface improves workflow, allowing clinicians to find controls quickly with fewer touches. Imaging is optimized for each exam type to minimize adjustments, improving time to diagnosis.

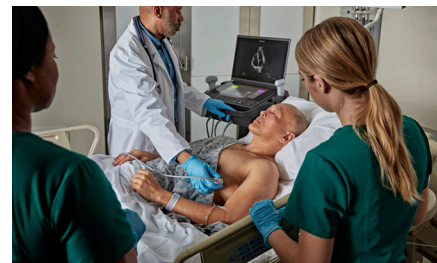

## Emergency Medicine

### Inclusive Workflow

Sonosite PX adapts to the way you provide emergency care. The streamlined user interface improves workflow, with touch and physical buttons providing access to the most-used controls for improved efficiency. The T8-3 transesophageal transducer paired with our cardiac resuscitation and lung exams increases capabilities for physicians working at the leading edge of ultrasound in Emergency Medicine.



## Musculoskeletal

### Confidently diagnose and treat patients in real-time.

Immediate in-office visualization helps physicians connect with patients to demystify treatment needs and improve patient satisfaction. Sonosite PX supports accurately visualizing fine anatomical details and assisting physicians to confidently diagnose and treat patients in real-time. With a single system interaction, the Auto Steep Needle Profiling helps physicians accurately guide needles for more precise injections, including injections into the delicate spaces between joints.



### Why Sonosite PX:

Sonosite ensures the security and privacy of its devices as well as the integrity of patients' Protected Health Information (PHI). The Manufacturer Disclosure Statement for Medical Device Security (MDS2) provides us a way to communicate with you in a standardized language for security. Whenever a new cyber security vulnerability is identified, we'll keep you informed of any potential impact to Sonosite systems through updates on our website. Sonosite PX meets FIPS199, FIPS200, and NIST 800-53.

**FUJIFILM SONOSITE**  
Value from Innovation

**FUJIFILM Sonosite, Inc.**  
Worldwide Headquarters  
21919 30th Drive SE, Bothell, WA 98021-3904  
Tel: +1 (425) 951 1200 or +1 (877) 657 8050  
Fax: +1 (425) 951 6800 [sonosite.com](http://sonosite.com)

**FUJIFILM Sonosite BV**  
European Headquarters  
Joop Geesinkweg 140, 1114 AB Amsterdam, Netherlands  
Tel: +31 20 462 0000

#### Sonosite Worldwide Offices

FUJIFILM Sonosite Australasia Pty Ltd: Australia . . . . . 1300 663 516  
FUJIFILM Sonosite Australasia Pty Ltd: New Zealand . . .0800 888 204  
FUJIFILM Sonosite Canada Inc. . . . . +1 888 554 5502  
FUJIFILM Sonosite GmbH-Germany . . . . . +49 69 80 88 40 30  
FUJIFILM Sonosite Iberica SL-Spain . . . . . +34 91 123 84 51  
FUJIFILM Sonosite India Pvt Ltd . . . . . +91 124 288 1100

FUJIFILM Sonosite Ltd-United Kingdom . . . . . +44 1462 341 151  
FUJIFILM Sonosite SARL-France . . . . . +33 182 880 702  
FUJIFILM Asia Pacific Pte. Ltd. . . . . +65 6383 9933  
FUJIFILM (China) Investment Co., Ltd . . . . . +86 21 5010 6000  
FUJIFILM Medical Co., Ltd.-Japan . . . . . +81 3 0418 7190  
FUJIFILM Middle East (FZE) . . . . . +971 4 887 8722

**Any patient. Anywhere. Anytime.**